




PRODUCT CATALOG // FIRE & PRESS // 2026

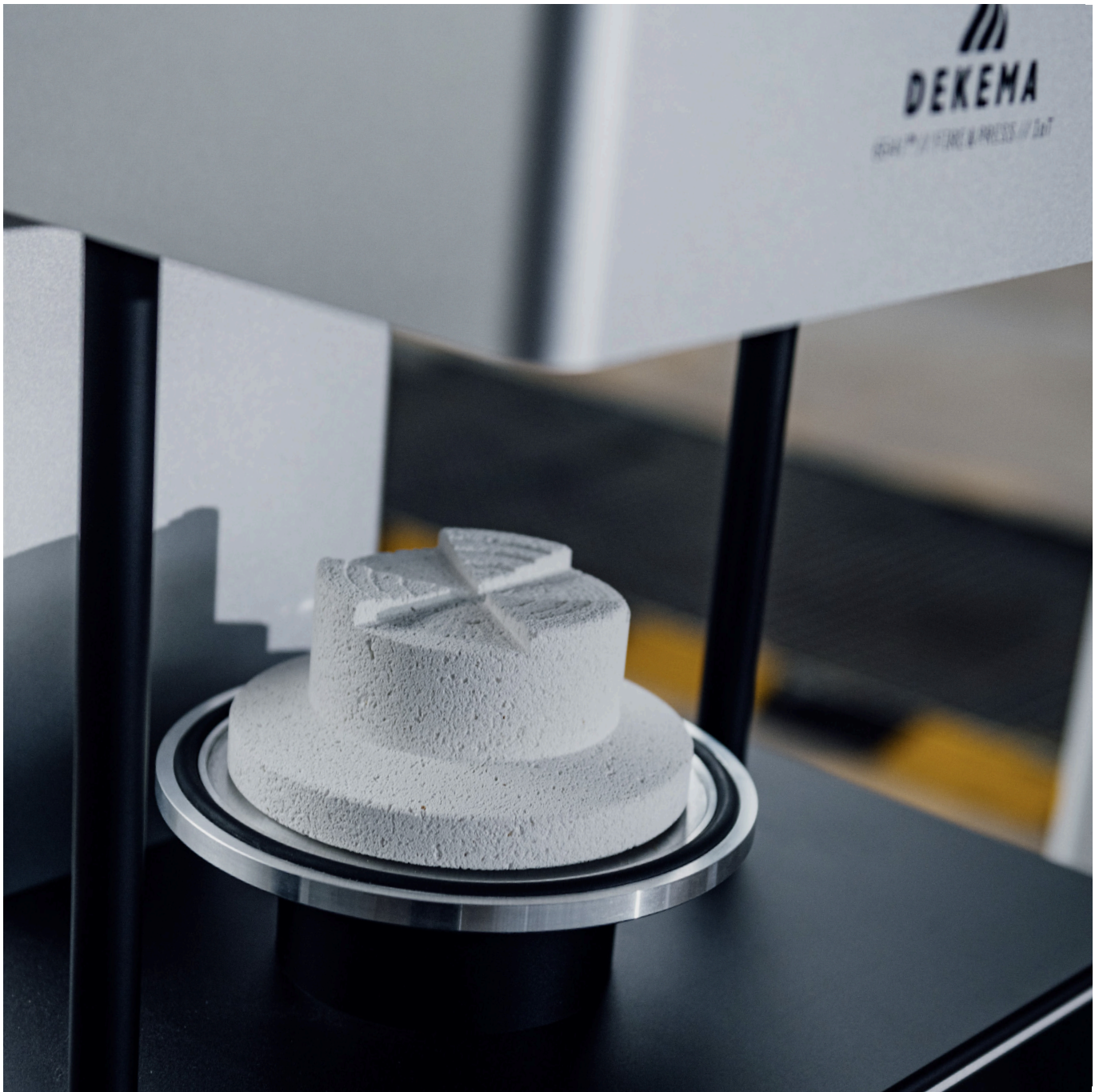

AUSTRIAN DPMAT™

FIRE + PRESS

654i™ G2

AUSTROMAT™ 654i™ G2

One for all – firing and pressing in perfection.




AUSTROMAT™ 654i™ G2

- » Even heat distribution through closed chamber design
- » Automatic control of the drying and cooling to optimize aesthetics
- » Multiple Press Time Automatic for optimal pressing results
- » Besides all common ring systems, press up to 20 g with the exclusive DEKEMA trixpress™ ring system
- » Press up to 3 different shades in one go, optimizing throughput
- » Even heat distribution over all surfaces of the investment ring which provides for exceptional aesthetics
- » DEKEMA AutoDry™
- » Access to more than 100.000 firing, pressing and sintering programs of all manufacturers of dental ceramics through the DEKEMA online

DEKEMA
Colorite™
You are the artist.

CREATE YOUR

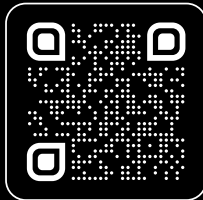

FURNACE

Applications	with vacuum	without vacuum
Fire	»	»
Press	»	»
Sinter	–	–
Glaze	»	»
Crystallize	»	»
Infiltrate	»	»



DEKEMA. Press.

Pressing shapes high-precision dentures. A blank is heated and pressed into the desired shape – for aesthetic and long-lasting results.

Unrivalled worldwide

The DEKEMA trixpress™ system combines an innovative upward pressing process with 3-plunger technology. This enables the processing of up to three colours in one step.

DEKEMA. Accessories.

AUSTROMAT™ 654i™ G2

Press table, multiple plunger press system Art. No. 101726

Press table for use in press furnaces with multiple plunger press system.



Firing table, multiple-plunger press furnaces Art. No. 101734

Firing table for use in press furnaces with multiple-plunger press system.



Ceramic supports (6 pieces) Art. No. 105488

Ceramic supports (6 pieces) for decoupled placement of press rings inside the preheating furnace.




CALIBRITE™ press set Art. No. 101818-01

Silver wire calibration set for exact temperature adjustment of press furnaces with multiple plunger press system.



CALIBRITE™ press wire Art. No. 101822-01

15 pre-formed silver wires for precise temperature calibration of press furnaces with multiple plunger press system.

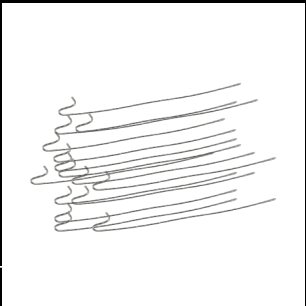

Firing tray, set Art. No. 106490

Firing tray, set of temperature-optimized honey-comb tray and 10 ceramic pins for exact positioning of firing objects.



DEKEMA. Accessories.

AUSTROMAT™ 654i™ G2

Tweezers, DEKEMA logo

Art. No. 105236

Furnace tweezers for safe handling of firing objects.



trixpress™ press ring, 380 g segments set

Art. No. 107210

Press ring for pressing over extensive restorations.



trixpress™ ring system, set

Art. No. 101807-01

Investment ring system trixpress™, set, consisting of base plate POM, segments for 100g and 200g investment rings, lid, 2 aluminum holding rings, Placeholder (50 pieces in blue color, Ø 13 mm) and 5 disposable press plungers.



Placeholder, blue color, Ø 13 mm (50 pieces) Art. No. 101915

Placeholder (50 pieces in blue color), Ø 13 mm, combustible without residue, for trixpress™ investment ring system, for sprewing press objects.



trixpress™ disposable plungers, Ø 13 mm (50 pieces)

Art. No. 105891

Disposable press plunger (50 pieces), Ø 13 mm, for single and multiple plunger pressings of all common investment materials and pressable ceramics.



DEKEMA. Technical data.

	AUSTROMAT™ 220	AUSTROMAT™ 224	AUSTROMAT™ 624i™ G2	AUSTROMAT™ 654i™ G2
Size (WxHxD)	10,47" x 20,91" x 10,16"	10,47" x 20,91" x 10,16"	15,08" x 25,71" x 13,50"	15,08" x 25,71" x 13,50"
Weight	28,5 lbs	28,5 lbs	42 lbs	53 lbs
Power supply	95 – 105 V/50–60 Hz or 110 – 125 V/50–60 Hz or 220 – 240 V/50–60 Hz	95 – 105 V/50–60 Hz or 110 – 125 V/50–60 Hz or 220 – 240 V/50–60 Hz	95 – 105 V/50–60 Hz or 110 – 125 V/50–60 Hz or 220 – 240 V/50–60 Hz	95 – 105 V/50–60 Hz or 110 – 125 V/50–60 Hz or 220 – 240 V/50–60 Hz
Maximum power supply	1.500 W	1.500 W	1.500 W	1.500 W
T _{max}	1.100 °C	1.100 °C	1.200 °C	1.200 °C

DEKEMA. Device equipment.

	AUSTROMAT™ 220	AUSTROMAT™ 224	AUSTROMAT™ 624i™ G2	AUSTROMAT™ 654i™ G2
Housing: Aluminum				
Corporate Design Customization	–	–	optional	optional
Touchscreen	4,3", color	7,0", color	10,4", color	10,4", color
USB interface/Ethernet interface	–/–	 /–		
Quartz glass heating element: DEKEMA specific heating element coil guided in pure quartz glass				
Fiber insulation: DEKEMA specific fiber insulation				–
DEKEMA stone insulation: internally fabricated, specific stone insulation	–	–	–	
Vacuum pump – version DEKEMA 838.3 R	optional	–		
Vacuum pump – DEKEMA AERIS™	–		–	–
DEKEMA ventilation unit	optional	optional	integrated	integrated
DEKEMA AutoDry™ & AUTORITE™	–	–		
MPTA (Multiple Press Time Automatic)	–	–	–	

DEKEMA. Software features.

	AUSTROMAT™ 220	AUSTROMAT™ 224	AUSTROMAT™ 624i™ G2	AUSTROMAT™ 654i™ G2
DEKEMA Internet database access for registered users				
DEKEMA Internet database furnace connectivity	–	–		
Internet program streaming	–	–		
DEKEMA Circle: Data Management in the Cloud	–	–		
IoT (Internet of Things) enables intelligent networking and real-time monitoring of the furnace – ensuring 99,9 % availability and seamless operation through preventive service.	–	–		
myDEKEMA: VNC-Viewer for remote control	–			
QM (quality management): Record of process files	–	–		
USB: Transferring firing programs (dds)	–			
USB: Transferring fotos (jpeg)	–	–		
USB: Storage of QM files (dfl)	–	–		
USB: Storage of Printouts	–	–		
USB: Backup / Update	–			
Internet: Backup / Update	–	–		
FTP-Server	–	–		
FTP-Client	–	–		
VNC-Server: VNC viewer for remote control of the DEKEMA furnace on the PC	–	–		
OPC: Industrial standardised connection for data exchange between different systems	–	–		
OPC UA: Secure, manufacturer-independent interface for Industry 4.0 and IoT communication	–	–	optional	optional
Programing languages	–	–	Classic, Professional	Classic, Professional, Press
Language-independent program management			–	–

DEKEMA. Accessories.

	AUSTROMAT™ 220	AUSTROMAT™ 224	AUSTROMAT™ 624i™ G2	AUSTROMAT™ 654i™ G2
Firing table, fibre (Art. No. 101733)				–
Firing tray, set (Art. No. 106490)	optional	optional	optional	optional
CALIBRITE™ fire set (Art. No. 101817-01)	optional			–
Firing table, multiple-plunger press furnaces (Art. No. 101734)	–	–	–	
Press table, multiple plunger press system (Art. No. 101726)	–	–	–	
Ceramic supports, 6 pieces (Art. No. 105488)	–	–	–	optional
CALIBRITE™ press set (Art. No. 101818-01)	–	–	–	
trixpress™ ring system, set (Art. No. 101807-01)	–	–	–	optional
trixpress™ press ring, 380 g segments set (Art. No. 107210)	–	–	–	optional
trixpress™ disposable plungers, Ø 13 mm, 50 pieces (Art. No. 105891)	–	–	–	optional
Placeholder, blue color, Ø 13 mm, 50 pieces (Art. No. 101915)	–	–	–	optional
Tweezers, DEKEMA logo (Art. No. 105236)	optional			
DEKEMA Interlock™ unit (Art. No. 107299)	–	optional	–	–
trix™4all, fully automated 3D nesting software for integration of third-party 3D printers into the trix™ workflow, 5-year license (Art. No. 107023)	–	–	–	optional



DEKEMA Dental-Keramiköfen GmbH
Industriestraße 22 • 83395 Freilassing • Germany
☎ +49-8654-46390 📠 +49-8654-66195
info@dekema.com • www.dekema.com

